

WEIGHT CAPACITY:

1000 lbs.

DIMENSIONS:

72" x 18" x 1.75"

WEIGHT:

18 lbs.

FEATURES:

100% X-Ray Translucent,
MRI and CT Scan compatible
250 lbs. Float Capacity

DESIGN:


1 Rotomolded one piece
seamless HDPE shell over non-toxic
Eco-Friendly polyurethane foam*
Fourteen over sized raise hand holds
Available with or without
12 or 22 speed clip pins

RELATED PRODUCTS:

Straps & Strapping Systems
Sof Loc Head Immobilizer

STANDARD COLORS:

Yellow, Blue, Orange, Red
Lime Green, Olive Drab, Black

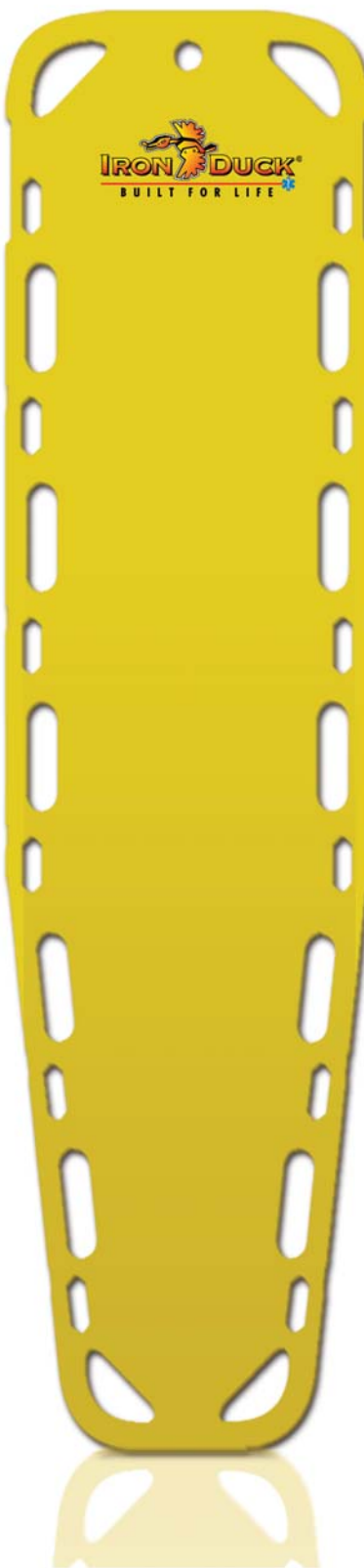

Ultra Vue 18

Product#: 35755

With 12 Pins: 35755-P

With 22 Pins: 35755-EP

Test results prove that Iron Duck's Ultra-Vue 18" backboard is among the most rigid boards in any weight class, outperforming competitors at weights up to 1,000 lbs.

In addition to the lowest deflection rates at any weight, the Vue 18 is X-ray translucent, MRI and CT scan compatible. Flat patient surface won't compromise CPR efforts or aggravate cervical spine injuries (as concave boards may).

Fabricated with rotomolded HDPE shell over nontoxic eco-friendly foam, the subtle surface texture prevents slippage but is easy to clean and disinfect. The Vue 18 ships with 10 two inch long reflective strips for added safety in low light situations.

*Kyoto & Montreal Compliant.

© Copyright 2013. Iron Duck, Built for Life, and Iron Duck logo are registered trademarks of Fleming Industries Inc. All Rights Reserved.